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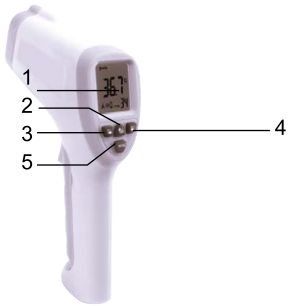
THERMOMETER

QUICK START GUIDE

1

OPIS

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GOCLEVER Sp. z o.o.

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EN**DESCRIPTION**

This product is a non-touch thermometer that measures both objects' and body temperature. The device is only for personal use. It does not substitute doctor's examination. The device stores up to 32 results and alarm for higher body temperature (over 38 °C).

DEVICE (PIC 1)

1. LCD display
2. backlit/laser key
3. SET key
4. ADJ key
5. MODE key
6. laser
7. transducer
8. on/measurement key
9. battery lid

TEMPERATURE MEASUREMENT

In order to measure body temperature, make sure that there is BODY on the display. In order to measure object's temperature, make sure that there is SURFACE on the display. You can change the BODY/SURFACE mode by pressing MODE key (Pic 1-5) as long as you get BODY or SURFACE.

YOU CAN MEASURE THE TEMPERATURE:

- » for body temperature: by aiming the laser light at a patient's forehead and pressing on/measurement key (Pic 1-8)
 - » for object temperature: by aiming the laser light at an object and pressing on/measurement key (Pic 1-8)
- The distance between measured object and a thermometer should be 5–15 cm.

SETTINGS

Press SET (Pic 1-3) for a few seconds to get into settings :

- » press one time – changing temperature unit °C/F
- » press two times – changing correction value (it corrects the result)
- » press three times – sound on/off

Change the options by pressing ADJ (Pic 1-4)

MEMORY

In order to see the results saved in device's memory, press ADJ (Pic 1-4) and SET (Pic 1-3) (up/down). In order to delete a particular result, press ADJ for a few seconds. You can remove all of the results from the memory by pressing ADJ and SET at the same time.

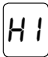
BATTERY INSTALLATION

1. Open battery lid (Pic 1-9) by pushing it.
2. Put white ribbon down.
3. Insert one AA battery inside (battery pole + up).
4. Insert second battery AA (battery pole + down)
5. Close battery lid.
6. Wait about 10 minutes before use after installing new batteries.

SAFETY PRECAUTIONS

1. Only use soft cotton cloth, dry or lightly damp with water or mild liquid, to clean the device.
2. Do not aim laser light at patient's eyes.
3. Keep the thermometer away from dampness and excessive sunlight.
4. Make sure before measurement of body temperature that there are no hair laying on the forehead and the skin is not covered with sweat.
5. Keep in mind that the results are only for reference for the preliminary diagnosis, you should always consult any worrying symptoms to your doctor.
6. Take away the batteries from the device if it has not been used for a long time
7. Only use standard batteries advised by a supplier.

TROUBLESHOOTING

Error	Cause	Solution
the display shows: 	body temperature is higher than 43 °C	try again
the display shows: 	body temperature is lower than 32 °C	try again

The errors above can be the result of the factors below:

- » body temperature has been measured incorrectly because of the contact with hair or sweat
- » too big distance between the thermometer and measured object/ patient's forehead
- » the measurements are being processed too soon one after another

SYMBOLS DEFINITION

	Attention! Please refer to this manual.
	Complies with MDD93/42/EEC requirements
	Type BF equipment (Refer to IEC 60601-1:1995)
	ENVIRONMENT PROTECTION – Waste electrical products should not be disposed of with household waste. Please recycle where facilities exist. Check with your local authority or retailer for recycling advice
	Keep dry

CORRECT DISPOSAL OF THIS PRODUCT

(Waste Electrical & Electronic Equipment)

(Applicable in the European Union and other European countries with separate collection systems)

This marking indicates that at the end of its working life this product should not be disposed of with other household waste in the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, please separate this from other types of wastes and recycle it responsibly to promote the sustainable reuse of material resources.

To dispose of your used device, please use the return and collection systems available in your area or contact the retailer where the product was purchased. They can take this product for safe environmental recycling.

SPECIFICATION

Product type	Thermometer
Measuring method	non-touch
Memory	32 results
Measuring modes	body temperature, surface temperature
Measuring range	32-43 °C (93.2-109 °F)
Accuracy	+/-3 °C (+/- 0.54 °F)
Measuring distance	5-15 cm
Power	2 AA size alkaline batteries (not included)

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CINDA INSTRUMENT CO. LTD, B3-913 WANDA PLAZA,